

# OptiPrep™ Application Sheet S50

## Fractionation of acidocalcisomes from trypanosomes (and other protozoa) and Gram-negative bacteria

- In Gram-negative bacteria these particles are often termed “volutin granules”
- The methodology can also resolve contractile vacuoles from protozoa (see Section 2c-2)
- ◆ OptiPrep™ is a 60% (w/v) solution of iodixanol in water, density = 1.32 g/ml
- ◆ To access other Application Sheets referred to in the text: return to the **2020SMemapp** file and select the appropriate S-number.
- ◆ Important technical notes, information regarding the extension of the technology to a variety of organisms and membrane analysis are contained in the “Technical Notes and Review” (Section 5).
- ◆ See also Application Sheet S58 for other organelle isolation methods from protozoa
- ◆ **IMPORTANT NOTE: OptiPrep™ Reference List RS14 gives an up-to-date bibliography of all the published papers reporting the use of iodixanol gradient methodology for the purification of organelles from non-mammalian eukaryotes. For yeast see Reference List RS15. To access these two files return to the initial list of Folders and select “Reference Lists”.**

### 1. Background

Acidocalcisomes, the electron-dense, acidic, calcium-storing organelles, which are rich in calcium and polyphosphate, were originally identified in *Trypanosoma cruzi* by Docampo et al [1] in 1995. They have now been detected in a range of prokaryotic and eukaryotic cells but there are significant differences in their characteristics from organism to organism. Although strict parallels in mammalian and non-mammalian cells may not exist, some dense granules in human platelets [2] and sea urchin eggs [3] share the same high levels of calcium and polyphosphate and may have similar functions.

Scott and Docampo [4] developed a discontinuous iodixanol gradient for purifying acidocalcisomes that effectively replaces the previous Percoll® method and overcomes many of the disadvantages associated with the use of Percoll®. The more general problem, observed with the recovery of mammalian organelles (e.g. lysosomes and peroxisomes) free from the colloidal silica particles with which they tend to co-sediment, is more serious in the case of acidocalcisomes because of their higher density and more rapid rate of sedimentation. Iodixanol, being a true solute, does not pose this problem at all. Moreover, since it is non-light-scattering, most spectrophotometric assays and electrophoresis can be carried out directly on gradient fractions without the need to remove the medium, so long as the concentration of the organelle in the gradient fraction is sufficiently high for accurate analysis. In some cases the gradient is able to provide simultaneous resolution of some other organelles.

### 2. Trypanosomes and other protozoa

#### 2a. Solutions required (see Section 2d-1)

- A. OptiPrep™ (60% w/v, iodixanol)
- B. Lysis buffer: 0.125 M sucrose, 50 mM KCl, 4 mM MgCl<sub>2</sub>, 0.5 mM EDTA, 5 mM dithiothreitol (DTT), 20 mM Hepes-KOH, pH 7.2
- C. OptiPrep™ diluent: 0.125 M sucrose, 0.3 M KCl, 24 mM MgCl<sub>2</sub>, 3.0 mM EDTA, 30 mM DTT, 120 mM Hepes-KOH, pH 7.2
- D. Iodixanol (50% w/v) working solution: mix 5 vol. of OptiPrep™ + 1 vol. of Solution C

Keep the following stock solutions at 4°C:

1 M Hepes (free acid)	23.8 g per 100 ml water
1 M KCl	7.45 g per 100 ml water
1 M MgCl <sub>2</sub> •6H <sub>2</sub> O	20.3 g per 100 ml water
100 mM EDTA(K <sub>2</sub> )	3.68 g per 100 ml water
100 mM DTT	1.54 g per 100 ml water

Solution B: Dissolve 8.55 g sucrose in 100 ml water; add 4 ml, 10 ml, 0.8 ml, 1 ml and 10 ml respectively of Hepes, KCl, MgCl<sub>2</sub>, EDTA and DTT stock solutions; adjust to pH 7.2 with 1 M KOH and make up to 200 ml.

Solution C: Dissolve 4.28 g sucrose in 30 ml of KCl stock solution; add 12 ml, 2.4 ml, 3 ml and 30 ml respectively of Hepes, MgCl<sub>2</sub>, EDTA and DTT stock solutions; adjust to pH 7.2 with 1 M KOH and make up to 100 ml.

**2b. Ultracentrifuge rotor requirements**

Ultracentrifuge with swinging-bucket rotors to accommodate approx 30 ml thick-walled tubes, e.g. Beckman SW28, Sorvall AH629 or similar (see [Section 2d-2](#))

**2c. Protocols**

Carry out all operations at 0-4°C.

**2c-1. Isolation of acidocalcisomes (adapted from ref 4)**

1. Lyse the washed epimastigotes by grinding with silicon carbide using standard procedures (see ref 3 for further details).
2. Centrifuge the lysed cells in Solution B at 144 g for 5 min (see [Section 2d-3](#)).
3. Decant the supernatant and centrifuge this at 325 g for 10 min (see [Section 2d-3](#)).
4. Decant and retain the supernatant.
5. Resuspend the pellet in Solution B and repeat the centrifugation at 325 g for 10 min.
6. Combine the two supernatants and centrifuge at 10,500 g for 30 min (see [Section 2d-4](#)).
7. Resuspend the pellet in Solution B (4 ml) by repeated passages through the 22-gauge needle of a syringe (see [Section 2d-4](#)).
8. Prepare density gradient solutions containing 24%, 28%, 34%, 37% and 40% (w/v) iodixanol by diluting Solution D with Solution B (see [Section 2d-5](#)).
9. In tubes for the chosen rotor, layer 4 ml of each of the density gradient solutions and layer the resuspended 10,500 g pellet on top of the discontinuous gradient, to fill the tube (see [Section 2d-6](#)).
10. Centrifuge at 50,000 g for 60 min.
11. The acidocalcisomes form a pellet at the bottom of the gradient. If the primary interest is isolation of the acidocalcisomes, aspirate the gradient and resuspend the pellet in a suitable medium. If the aim is to analyze other organelles as well, then collect the gradient in a series of 15-20 equal volume fractions (see [Sections 2d-7 and 2d-8](#)).

**2c-2. Isolation of acidocalcisomes and contractile vacuoles (adapted from ref 5)**

1. Clarify the lysate sequentially using two centrifugations at 36 g for 5 min and one at 144 g for 10 min (see [Section 2d-3](#)).
2. Centrifuge the final supernatant at 100,000 g for 1 h.
3. Resuspend the pellet in Solution B (2 ml) by repeated passages through the 22-gauge needle of a syringe (see [Section 2d-4](#)).
4. Mix the suspension with an equal volume of Solution D and 0.15 ml of Solution B to adjust it to 24% (w/v) iodixanol.
5. Prepare density gradient solutions containing 15%, 20%, 28%, 34%, 37% and 40% (w/v) iodixanol by diluting Solution D with Solution B.

6. In tubes for the chosen rotor, layer 4 ml of each of the density gradient solutions, including the sample in 24% iodixanol ([see Section 2d-6](#)).
7. Centrifuge at 50,000 g for 60 min.
8. Collect the gradient in a series of 15-20 equal volume fractions ([see Section 2d-7 and 2d-8](#)).

## 2d. Technical Notes and Review

### 2d-1. Homogenization media

Solutions are commonly buffered with Tris, Hepes, Tricine or triethanolamine (at 10-20 mM concentration) and it is unlikely if the type of buffer significantly influences the fractionation.

The preparation of a Working Solution as described, ensures that the concentrations of KCl, MgCl<sub>2</sub>, EDTA, DTT and the buffer (Hepes-KOH, pH 7.2) are constant throughout the gradient. If this is deemed unimportant the iodixanol solutions may be prepared simply by diluting OptiPrep™ with Solution B. Strategies for preparing working solutions are given in [Application Sheet S01](#).

Protease inhibitors may be included in Solutions B and C at the operator's discretion.

### 2d-2 Ultracentrifuge rotors

The gradient + sample volume is 20-24 ml, thus thick-walled tubes which may be partially filled are the recommended ones for the rotor (total tube volume = approx. 30 ml). There is no obvious reason however why the gradient cannot be scaled up to allow the use of thin-walled tubes, which may be more convenient if complete unloading of the gradient into multiple fractions is envisaged. Similarly the procedure may be scaled down for use in rotors with a smaller tube volume.

### 2d-3 Clarification of the lysate

In many cases of acidocalcisome preparation, for example from *Toxoplasma gondii* [6-8], the 144 g and 325 g steps have been replaced by 36 g and 144 g.

### 2f-4 Preparation of a crude acidocalcisome fraction

In the case of *Toxoplasma gondii* [6,7] the crude acidocalcisomes were sedimented at 15,000 g rather than 10,000 g. Suspension of the crude fraction should be carried out as gently as possible to avoid damage not only to the organelles of interest but also to any other organelles present – particularly those which may release degradative enzymes. If median loading of the sample in the gradient is chosen, rather than top-loading, suspend the pellet in no more than approx 2 ml of Solution B so the volume after adjustment of the density with Solution D remains manageable.

### 2d-5 Gradient format (all iodixanol concentrations are % (w/v))

In some instances, for example in the isolation of acidocalcisomes from *Dictyostelium discoideum* the 40% iodixanol layer has been omitted [9]. To avoid any possible loss of material due to the rapid accumulation of particles at the sample/24% iodixanol interface in [Protocol 2c-1](#), it may be preferable to suspend the crude pellet in one of the gradient layers, e.g. 24% or 28% iodixanol as in Protocol 4b. In [Protocol 2c-2](#) the lower density layers (15% and 20% iodixanol) are included to improve the resolution of the contractile vacuole from the denser mitochondria, glycosomes and lysosomes; such a format may generally be beneficial in any studies in which a more complete fractionation is required.

For *Toxoplasma gondii* the gradient comprised 10%, 15%, 20% (sample), 25% and 30% iodixanol [6,7], while for isolation of acidocalcisomes from the plant trypanosomatid *Phytomonas françai* [Protocol 2c-1](#) was used [10]. In proteomic analytical studies of acidocalcisomes, Ferella et al [11] used a discontinuous 20-50% (w/v) iodixanol gradient and in a similar study of the contractile vacuole the discontinuous gradient comprised 15, 20, 25, 30, 34, 37 and 40% iodixanol, with the sample in the 25% layer [12].

## 2d-6 Forming the discontinuous gradient

Although overlaying (i.e. starting with the densest layer) is the most common means of creating a discontinuous gradient, underlayering (i.e. starting with the least dense layer) with a syringe and metal cannula is more reliable and the recommended method for making discontinuous gradients. For more information on gradient construction see [Application Sheet S03](#). If necessary, adjust all volumes proportionately so that tubes (after sample application) are properly filled according to the manufacturer's instructions.


## 2d-7 Fractionating the gradient

If thick-walled tubes are used then aspiration from the bottom of the tube or from the meniscus (with a Labconco Auto Densi-flow device) are acceptable methods for harvesting gradients. Unloading by upward displacement with a dense liquid may be less acceptable since the pellet may become dispersed. So long as the acidocalcisomes do not form a too firmly-packed pellet, tube puncture may be satisfactory for a thin-walled tube. For more information on harvesting gradients see [Application Sheet S52](#).

## 2d-8 Analysis

Protocol 2c-1 has been used preparatively for the isolation of acidocalcisomes from *Trypanosoma cruzi* [4,13-15], *Trypanosoma brucei* [15] and *Leishmania major* [15] in studies on  $\text{Ca}^{2+}$  and phosphate metabolism. This protocol also permits the partial resolution of some other organelles [16], although glycosomes, lysosomes and the vacuolar compartment tend to overlap close to the top interface between the sample and the top layer of the gradient.

There are however quite clear differences in the manner in which organelles from different organisms behave using Protocol 2c-1. Figure 1 shows the distribution of organelles from a 10,500g pellet prepared from a post-nuclear supernatant of a *Dictyostelium discoideum* homogenate, which was layered on top of the gradient, as in Protocol 2c-1, but omitting the 40% iodixanol layer [9]. Alkaline phosphodiesterase shows a clear concentration at the top of the gradient (Figure 1); this is an established marker for contractile vacuoles from a number of microorganisms. The vacuole invariably bands at a lower density than any other organelle. On the other hand both pyrophosphatase and the vacuolar  $\text{H}^+$ -ATPase (V- $\text{H}^+$  ATPase) are present in both the vacuole and the dense acidocalcisome, confirming a functional link between these two particles [9]. Acid phosphatase (lysosomes) and succinate cytochrome c reductase (mitochondria) each show distinctive profiles. The distribution of the mitochondria is broad (compared to the usual pattern from mammalian cells); nevertheless fractions 8-10 which demonstrate the highest succinate-cytochrome c reductase are also significantly impoverished in markers for other organelles. It is also notable that the lysosomes from *Dictyostelium* are denser than are the bulk of the mitochondria (the reverse is true in all mammalian cells so far studied).

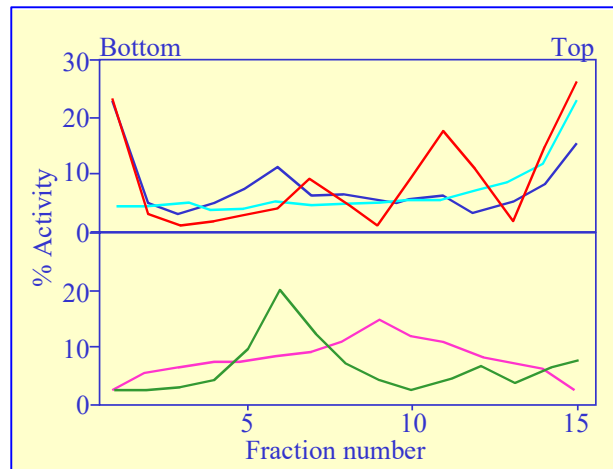

**Figure 1.** Fractionation of *Dictyostelium* 10,000g pellet: pyrophosphatase (—); V- $\text{H}^+$  ATPase (—), alkaline phosphodiesterase (—), acid phosphatase (—), succinate-cyt c reductase (—). Adapted from ref 8 with kind permission of the authors and The American Society of Biochemistry and Molecular Biology.

- ◆ In the case of *Chlamydomonas reinhardtii* Protocol 2c-1 was also able to resolve very clearly an acidocalcisome fraction and in this case the mitochondria were also well resolved about a third of the way down the gradient, but the chloroplasts were rather broadly distributed [17].

Protocol 2c-2 provides a much more clear resolution in the case of *Toxoplasma gondii* [6,7] and *Trypanosoma cruzi* [5] of the vacuole, lysosomes and glycosomes and the acidocalcisomes which band towards the top, middle and bottom of the gradient respectively.

- ◆ Other papers using this methodology for *Trypanosoma brucei* report that modulation of polyphosphate alters acidocalcisome biogenesis and function [18]; acidocalcisome lipids and glycolipids [19] and that the inositol 1,4,5-trisphosphate receptor present on acidocalcisomes has an important role in growth and infectivity [20].
- ◆ A discontinuous iodixanol gradient has also been used to purify contractile vacuoles from *Dictyostelium* [21].

### 3. Bacteria

#### 3a Clarification of the lysate

The lysates from *Rhodospirillum rubrum* and *Agrobacterium tumerfaciens*) were clarified at 1000g for 5 min [22,23].

#### 3b Preparation of a crude acidocalcisome fraction

For bacteria 14,500g was used to produce a crude organelle fraction [22,23]. Thus to maximize yields and purity of the acidocalcisomes it may be necessary to optimize the differential centrifugation of the homogenate.

#### 3c Gradient format (all iodixanol concentrations are % (w/v))

Bacteria suspensions (*Rhodospirillum rubrum* and *Agrobacterium tumerfaciens*) were adjusted to 30% iodixanol and made part of a 24%, 28%, 30%, 35%, 40% iodixanol gradient [22,23] Note that in some bacteria the acidocalcisomes may be termed volutin granules [23]. A 14,500 g (10 min) fraction from a *Corynebacterium matruchotii* lysate was also median loaded as the 30% (w/v) step in a 24%, 28%, 30%, 35%, 40% (w/v) iodixanol gradient and centrifuged at 235,000  $g_{av}$  for 1 h. The particles banded at the 30%-35% interface [24].

- ◆ A review of the functional significance of acidocalcisomes is to be found in ref 25
- ◆ A very thorough proteomic analysis of acidocalcisomes from a variety of sources, using variants of the above technology, see ref 26.

### 4. References

1. Docampo, R., Scott, D. A., Vercesi, A. E. and Moreno, S. N. J. (1995) *Intracellular Ca<sup>2+</sup> storage in acidocalcisomes of Trypanosoma cruzi* Biochem J., **310**, 1005-1012
2. Docampo, R., de Souza, W., Miranda, K., Rohloff, P. and Moreno, S.N.J. (2005) *Acidocalcisomes – conserved from bacteria to man* Nature Rev. Microbiol., **3**, 251-261
3. Ramos, I.B., Miranda, K., Pace, D.A., Verbist, K.C., Lin, F-Y., Zhang, Y., Oldfield, E., Machado, E.A., de Souza, W. and Docampo, R. (2010) *Calcium- and polyphosphate-containing acidic granules of sea urchin eggs are similar to acidocalcisomes, but are not the targets for NAADP* Biochem. J., **429**, 485–495
4. Scott, D. A. and Docampo, R. (2000) *Characterization of isolated acidocalcisomes of Trypanosoma cruzi* J. Biol. Chem., **275**, 24215-24221
5. Rohloff, P., Montalvetti, A. and Docampo, R. (2003) *Acidocalcisomes and the contractile vacuole complex are involved in osmoregulation in Trypanosoma cruzi* J. Biol. Chem., **279**, 52270-52281
6. Rodrigues, C.O., Ruiz, F.A., Rohloff, P., Scott, D.A. and Moreno, S.N.J. (2002) *Characterization of isolated acidocalcisomes from Toxoplasma gondii Tachyzoites reveals a novel pool of hydrolysable polyphosphate* J. Biol. Chem., **277**, 48650-48656
7. Rohloff, P., Miranda, K., Rodrigues, J.C.F., Fang, J., Galizzi, M., Plattner, H., Hentschel, J. and Moreno, S.N.J. (2011) *Calcium uptake and proton transport by acidocalcisomes of Toxoplasma gondii* PloS One **6**: e18390
8. Ferreira, D. da S., Menezes Resende, I.T. and Lopez, J.A. (2014) *Proteome investigation of an organellar fraction of Toxoplasma gondii: a preliminary study* BMC Proc., **8 (Suppl 4)**: P74
9. Marchesini, N., Ruiz, F.A. Vieira, M. and Docampo, R. (2002) *Acidocalcisomes are functionally linked to the contractile vacuole of Dictyostelium discoideum* J. Biol. Chem., **277**, 8146-8153

10. Miranda, K., Rodrigues, C.O., Hentchel, J., Vercesi, A., Plattner, H., de Souza, W. and Docampo, R. (2004) *Acidocalcisomes of Phytomonas francai possess distinct morphological characteristics and contain iron* Microsc. Microanal., **10**, 647-655
11. Ferella, M., Nilsson, D., Darban, H., Rodrigues, C., Bontempi, E.J., Docampo, R. and Andersson, B. (2008) *Proteomics in Trypanosoma cruzi - localization of novel proteins to various organelles* Proteomics, **8**, 2735-2749
12. Ulrich, P.N., Jimenez, V., Park, M., Martins, V.P., Atwood III, J., Moles, K., Collins, D., Rohloff, P., Tarleton, R., Moreno, S.N.J., Orlando, R. and Docampo, R. (2011) *Identification of contractile vacuole proteins in Trypanosoma cruzi* PLoS One **6**: e18013
13. Martinez, R., Wang, Y., Benaim, G., Benchimol, M., de Souza, W., Scott, D. A. and Docampo, R. (2002) *A proton pumping pyrophosphatase in the Golgi apparatus and plasma membrane vesicles of Trypanosoma cruzi* Mol. Biochem. Parasitol., **120**, 205-213
14. Rohloff, P., Rodrigues, C. O. and Docampo, R., (2003) *Regulatory volume decrease in Trypanosoma cruzi involves amoni acid efflux and changes in intracellular calcium* Mol. Biochem. Parasitol., **126**, 219-230
15. Moreno, B., Urbina, J.A., Oldfield, E., Bailey, B.N. Rodrigues, C.O. and Docampo, R. (2000) *<sup>31</sup>P NMR spectroscopy of Trypanosoma brucei, Trypanosoma cruzi, and Leishmania major* J. Biol. Chem., **275**, 28356-28362
16. Ruiz, F.A., Rodrigues, C.O. and Docampo, R. (2001) *Rapid changes in polyphosphate content within acidocalcisomes in response to cell growth, differentiation and environmental stress in Trypanosoma cruzi* J. Biol. Chem., **276**, 26114-26121
17. Ruiz, F.A., Marchesini, N., Seufferheld, M., Govindjee and Docampo, R. (2001) *The polyphosphate bodies of Chlamydomonas reinhardtii possess a proton pumping pyrophosphatase and are similar to acidocalcisomes* J. Biol. Chem., **276**, 46196-46203
18. Fang, J., Rohloff, P., Miranda, K. and Docampo, R. (2007) *Ablation of a small transmembrane protein of Trypanosoma brucei (TbVTC1) involved in the synthesis of polyphosphate alters acidocalcisome biogenesis and function, and leads to a cytokinesis defect* Biochem. J., **407**, 161-170
19. Salto, M.L., Kuhlenschmidt, T., Kuhlenschmidt, M., de Lederkremer, R.M. and Docampo, R. (2008) *Phospholipid and glycolipid composition of acidocalcisomes of Trypanosoma cruzi* Mol. Biochem. Parasitol., **158** 120-130
20. Huang, G., Bartlett, P.J., Thomas, A.P., Moreno, S.N.J. and Docampo, R. (2013) *Acidocalcisomes of Trypanosoma brucei have an inositol 1,4,5-trisphosphate receptor that is required for growth and infectivity* Proc. Natl. Acad. Sci. USA, **110**, 1887-1892
21. Sivaramakrishnan, V. and Fountain, S.J. (2012) *A mechanism of intracellular P2X receptor activation* J. Biol. Chem., **287**, 28315-28326
22. Suefferheld, M., Lea, C.R., Vieira, M., Oldfield, E. and Docampo, R. (2004) *The H<sup>+</sup> -pyrophosphatase of Rhodospirillum rubrum is predominantly located in polyphosphate-rich acidocalcisomes* J. Biol. Chem., **279**, 51193-51202
23. Seufferheld, M., Vieira, M. C. F., Ruiz, F.A., Rodrigues, C.O., Moreno, S.N.J. and Docampo, R. (2003) *Identification of organelles in bacteria similar to acidocalcisomes of unicellular eukaryotes* J. Biol. Chem., **278**, 29971-29978
24. Linton, K.M., Tapping, C.R., Adams, D.G., Carter, D.H., Shore, R.C. and Aaron, J.E. (2013) *A silicon cell cycle in a bacterial model of calcium phosphate mineralogenesis* Micron, **44**, 419-432
25. Patel, S. and Docampo, R. (2010) *Acidic calcium stores open for business: expanding the potential for intracellular Ca<sup>2+</sup> signaling* Trends Cell Biol., **20**, 277-286
26. Huang, G., Ulrich, P.N., Storey, M., Johnson, D., Tischer, J., Tovar, J.A., Moreno, S.N.J., Orlando, R. and Docampo, R. (2014) *Proteomic analysis of the acidocalcisome, an organelle conserved from bacteria to human cells* PLoS Pathog., **10**: e1004555

OptiPrep™ Application Sheet S50: 10<sup>th</sup> edition, January 2020